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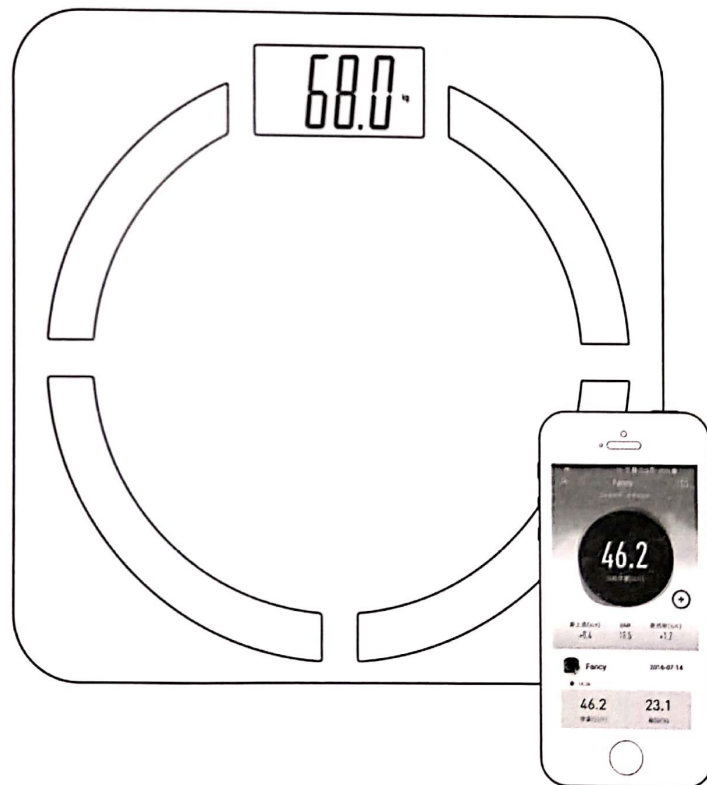
User Manual

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## Contents

### Instruction

Product instruction .....	1
Safety warning .....	1
Main parts .....	1
Battery installation .....	2

### Measuring

APP Installation .....	3
Registration .....	4
measuring .....	4
Measure data and report .....	5
Family members adding and recognition .....	5

### Features .....

	6
--	---

### recycling .....

	6
--	---

### Maintain .....

	7
--	---

### Warranty Claim .....

	8
--	---

### Contact Us .....

	9
--	---

## Instruction


### Product instruction:

Bluetooth scale helps you to store and keep track of your weight on mobile devices by OKOK app. This scale can transmit data of your Body fat, body hydration, body muscle, body bone, metabolism, BMI, Visceral Fat, Protein, Obesity Degree, Body Age, Weight without fat etc. to your smart devices via bluetooth.

### Safety warning :


1. Those who use heart pacemaker or with other implanted medical device should not use the device.  
That's because Bio-electric current through the body when measuring. Then the medical devices may not work normally.
2. The Fat%, hydration% content evaluated by the scale is only for reference (not for medical purpose).  
If your Fat%, hydration% content is over/below normal level, please consult your doctor for more advice.
3. People with below conditions are not advisable using the device; children under 10, elder above 80, and pregnant women.

### Main parts:
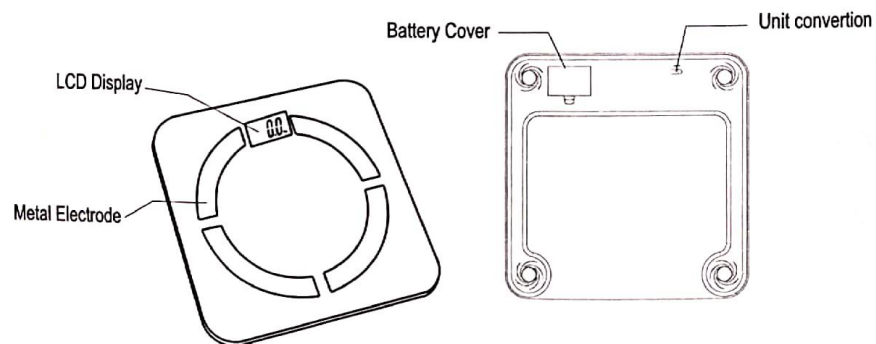
### Battery installation

1. Open the battery cover on the backside of the product.
2. Install the batteries and ensure + / - terminals are correct.
3. Close the battery door.

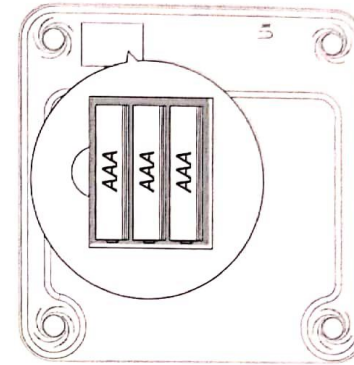Under below condition, please change new batteries

- Display **Lo**
- Dimmer display
- Nothing display

### ※ Note:

- Please use corresponding battery.
- Please remove battery if it is not being used for a long time.

## APP Installation

1. Compatible with the following mobile oSes:

IOS 7.0

Android 4.3


Bluetooth 4.0

2. Turn on your mobile phone's bluetooth and WIFI or 3G/4G.


3. Scan the QR below to install OKOK app.

also you can download OKOK app from the App store( for IOS only).
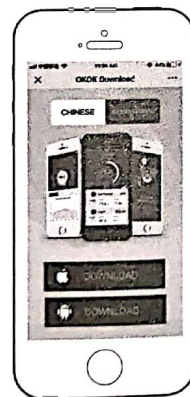
## Registration

Login

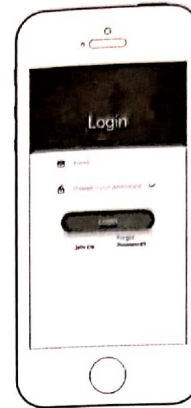profile setting

Target weight setting

measuring

Tips:


1. Step on the scale with bear feet.

Place the scale on a hard, flat platform

## Health data analyze & health tips

In the health score page, tap on every measurement, the analysis and health tip will come out

Body fat, body hydration, body muscle, body bone, metabolism, BMI, Visceral Fat, Protein, Obesity Degree, Body Age, Weight without fat
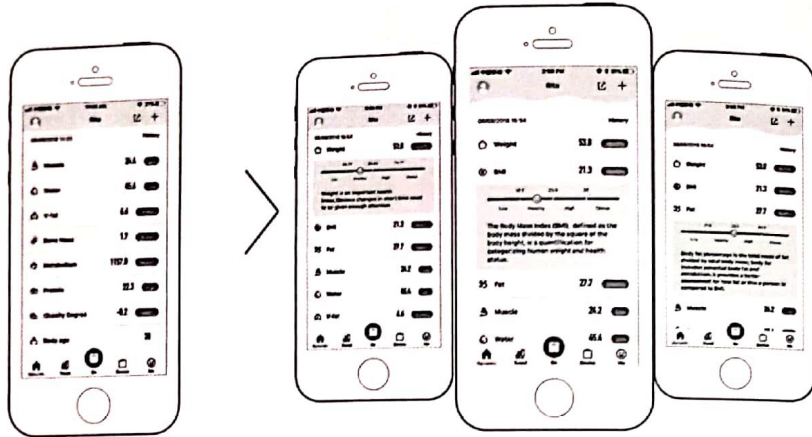
## Family members adding and recognition

Press "Me" on main page, to add more family members

A total of 8 family members can be added. Using a 3-kg difference in weight, the scale automatically recognise the family member.

## Features

Item	Parameter
Capacity	180kg/397lb/28st.5lb
Division	0.1KG/0.2lb/1lb
Unit	kg; lb; st:lb
High fat indication	ErrH
Low fat indication	ErrL
Height range	100 ~ 240/cm; 3' 3 ~ 7' 10/ft inch
Measuring body fat range	5% ~ 50%
Measuring water range	10% ~ 80%
Measuring bone range	0.5kg ~ 8.0kg
Measuring muscle range	10% ~ 80%
Auto identification range	3kg


## recycling

When no longer functional an electronic product should be disposed of with the least possible environmental impact. The devise should be disposed of according to your local municipal regulations, but in most cased you can take the product to your local recycling centre.

## Product maintain


Please follow below instruction to maintain the product, so that to keep it in good condition.

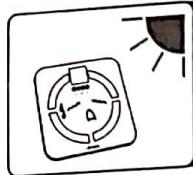Avoid to use the scale in humid, overheat or overcold environment.  
(temperature range from  $-20\sim 60^{\circ}\text{C}$ )

Do not put the product in water or wash it with water, to clean the scale with a slightly damp cloth when it is dirty.

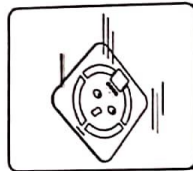Do not drop, shock or strike the scale as it is a precise instrument.

## Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from ITSU .
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non residential users.
8. All services covered by this warranty must be approved by ITSU and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, ITSU reserves the right to make substitution in lieu of repair or replacement.

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