ITSU

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User Manual
IS0141 使用手册

User Manual

All information is correct at the time of print

Thank you for purchasing ITSU Knee massager.

Please read the Safety Precautions carefully before using the foot massager. The Safety Precautions are divided into 'warnings' and 'attention' in terms of the degree of danger. This book should be kept available for future reference after reading.

CONTENTS

Safety Instructions	03
Warranty Claims	06
Product Parts	07
Product Features And Operation Illustrations	08
Product Operation Guide	09
Product Maintenance	10
Troubleshooting	11
Product Specifications	12
Contact Us	13

Safety Instructions

a) Important safety instructions

- To ensure a safe and proper operation, please read the instruction manual carefully before using this product.
- Use this product only as described in this manual.
- Always unplug this product from the electrical outlet immediately after use and before cleaning.
- Do not use this product if you have sensitive skin, neuropathy or poor circulation. Don't let children or people with mental disabilities use the heating function without proper supervision. Improper use may be dangerous.
- Do not operate this product under a blanket or when covered. Over-heating can occur and cause fire, electric shock or injury.
- Do not cover this product in order to avoid electrical faults or damages due to excessively high motor temperature.
- Stop the massage immediately if you feel any discomfort.
- Never operate this product if it has been dropped, has a damaged cord or plug; experiences a malfunction, or if it has been dropped in water. Please contact an authorized service centre to do the necessary check and repair.
- Do not use the components that are not recommended by ITSU.
- Do not insert any objects into any openings in the product.
- If you want to cut off the power supply, please make sure you have turned off all the operating controls, as well as remove the power plug from the socket.
- Please do not use this product as a medical device. If you need a professional diagnosis or treatment, please consult a professional medical person.

b) Operating Environment Precautions

- Do not use this massager in humid or extremely dusty environments to avoid malfunction.
- Do not use this massager at a temperature above 40 degree centigrade.
- Do not use this massager outdoors.
- Avoid placing or operating this massager in direct sunlight or near a heating appliance, stove.

Safety Instructions

c) User precautions

- Make sure that this product is far away from any hot surfaces.
- Please consult your doctor before use for those who have the following diseases:
 - Those who have heart disease, malignant tumor, back pain, spinal scoliosis, hypertension, osteoporosis, etc.
 - Those who are pregnant or having a menstrual period.
 - Those who have a physical discomfort, fever, a swelling or inflammation, bruises, allergies, etc.
 - Those who are warned by the doctor not to do any massage to avoid injury or physical discomforts such as blood clots (thrombosis), severe aneurysm, acute venous tumor, all kinds of skin inflammation, infectious skin disease (including subcutaneous tissue inflammation), etc.
 - Warnings: Do not use the massager together with the following medical instruments.
 - 1) Those who are using a pacemaker or other implantable medical equipment.
 - 2) Those who are using an artificial heart.
 - 3) Those who are wearing an electrocardiograph for similar medical equipment.

Otherwise it will cause malfunction for those medical equipment, putting the user's life in danger.

d) Precautions for power supply

- Do not turn on the switch when your hands are wet.
- Do not immerse in or spray water or liquid on the control power switch.

e) Adapter / USB cable Precautions

- Please match a suitable adapter by yourself, and then connect the Type-C cable that comes with this product.
- When charging, the charging indicator is red.
After being fully charged, the charging indicator will be green.
When the battery is low, the charging indicator flashes red with sound.

Safety Instructions

- Do not wind the power plug around the product lest it should be damaged and cause an electric shock or fire.
- Do not operate this product if it has a damaged cord or plug.
- Do not use this product if the power socket is loose.
- Do not carry this product by its supply cord or use the cord as a handle.
- When the power plug is damaged, please call the service centre or an authorized technical person for the necessary check and maintenance.




Attention

- Do not fall asleep while massager is on.
- Do not use this product while you are drunk or physically uncomfortable.
- Don't stand on the product or put a heavy weight on it.
- Do not use this product for more than 30 minutes.
- Do not use this product if the cloth cover is worn out.
- Do not use this product within an hour before or after a meal.
- When you are using the product at a low room temperature, do not suddenly use the heat function. You can adjust the room temperature to normal room temperature, before starting the heat function.
- When moving this product from a low temperature storage place to a warmer setting, it may cause malfunction due to water vapor gathered in the components. Set the product in a normal room temperature for one hour before using it.
- Do not use this product continuously if you find any components damaged or the product in an abnormal condition.
- When current interruption occurs during your massage, try to avoid further damage due to sudden power restoration.
- Stop using the product if there is thunder and lightning.

Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from ITSU.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by ITSU and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, ITSU reserves the right to make substitution in lieu of repair or replacement.

Product Parts

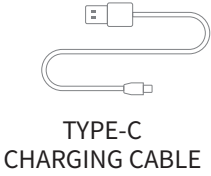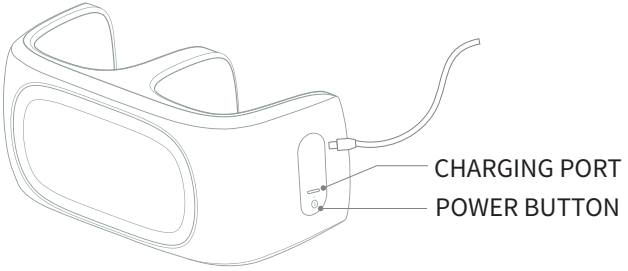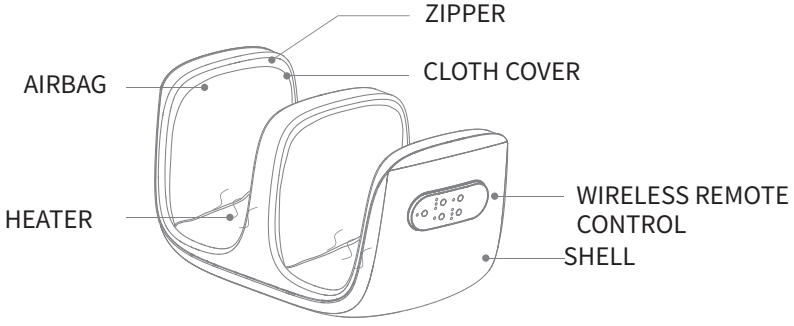
Product Features And Operaton Illustrations

Preparation before use

Thank you for purchasing this product. Before use, please read this manual carefully and pay special attention to safety matters to ensure correct operation and use. Please keep this manual properly for future reference.

Please place the product horizontally.


Instructions

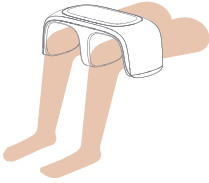LEG mode

KNEE mode

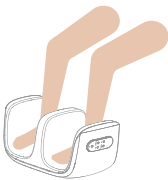FOOT mode

Product Operation Guide

Please switch on POWER button (on the product, located below remote controller) before use.

Function of remote controller's buttons:

ON/OFF button:

Press to turn on the massager

Press again to turn off .

(Default timer is 15 minutes)

Press to turn ON Auto 1 (flashing one).

Press again for Auto 2 (flashing twice).

Press again for Auto 3 (triple flashing).

Press to enable Leg mode.

Press again to enable Knee mode.

Press again to enable Foot mode.

press to switch airbag massage intensity
low-mid-high).

Press to enable heater.

Press again to disable heater.

- If a power failure occurs when the product is in use, turn off the power switch .
- Do not put hands or fingers close to and into the gap to avoid injury.
- Do not let children and pets play around the product to avoid accidents.
- Do not use the product at no load to avoid wearing the cloth cover.

Product Maintenance

a) Storage

- Please store it in the packing carton, and put in a safe, dry, well-ventilated, and shaded place.
- Do not twist the power plug.
- Do not hang the power wire and pull the product.
- Avoid sharp items from destroying the surface of the product.
- Please don't store the machine under direct sunlight or in places of high temperature.

b) Clean

- Before cleaning, make sure that the power cord is unplugged from the electrical outlet and that the machine is completely cool.
- Use a damp cloth to clean the product. Ensure the product is completely dry before reuse.
- The detachable fabric cloth can be washed in a normal washing way.

***Attention**

- Make sure not to use benzene or thinner and other corrosive detergents to clean the product.

c) Maintenance

- Keep the product clean and store in a dry environment.
- When the product is not in use for a long time, please put a dustproof cover on the product.

Troubleshooting

If you meet the following problems, please referring to the following solutions. If the product is still malfunctioning or doesn't work, please send the product to our service centre or an authorized technician for the necessary checks and repair.

Problem	Check points	Troubleshooting
Friction, air exhausting or noise from the compressor when the product is in use.	The noises occur from internal parts due to the mechanical design when the product is in use. They will not have any effect on the operation of the unit.	No need for any actions.
The unit does not start up.	Does the battery need to be charged? Check whether the POWER ON/OFF button has been turned on. Check whether no massage function has been selected. Protection from overheat.	Fully charge the battery before use. Turn on the power switch. Select a massage function. Turn off the power. Restart the unit after a 30 minute break.
Sudden cease of operation.	The rated operating time has elapsed (15 minutes).	Turn off the power. Restart the unit after a 10 minute break.

The noises occur from internal parts due to the mechanical design when the product is in use. They will not have any effect on the operation of the unit.

Product Specifications

Power supply	AC 220V-240V, 50/60Hz
Power consumption	12W
Timer	15 minutes
Dimension(LxWxH)	425 x 190 x 206 mm
Weight	2.6Kg
Dimension of packing box (LxWxH)	485 x 228 x 257 mm
Weight (including the color box)	3.1Kg

Contact Us

ITSU World Sdn Bhd
13 Lorong Teknologi B. Nouvelle Industry Park.
Taman Sains SG1 Kota Damansara. PJU 5.
47810 Petaling Jaya
Malaysia
Tel: +603-6732 3722

ITSU World(HK) Limited
Workshop D1,28/F. TML Tower
No. 3 Hoi Shing Road, Tsuen Wan, N.T.
Tel: +852 3611 2228

PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel: +62 21-6983 1123

User Manual

All information is correct at the time of print

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